



BON CHARGE



BASE STAND

USER GUIDE


1. COMPATIBILITY

This stand is compatible with Max Red Light Therapy Device and Super Max Red Light Therapy Device.


2. ASSEMBLY INSTRUCTIONS

Items included:

x2 BRACKETS

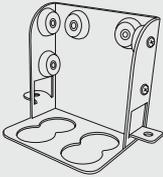


x1 BASE

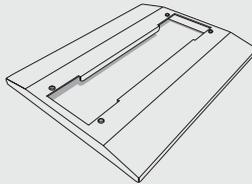


Step 1:

Unscrew the 4x bolts on the surface of the base stand.

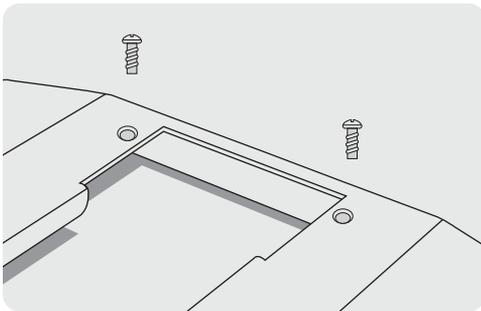


Step 2:

Place the brackets on the top of the screw holes.

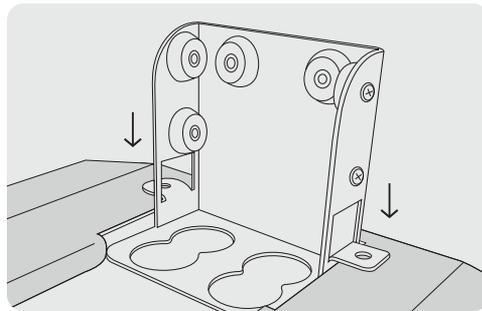


Step 3:

Screw the bolts back in place, repeat other side.

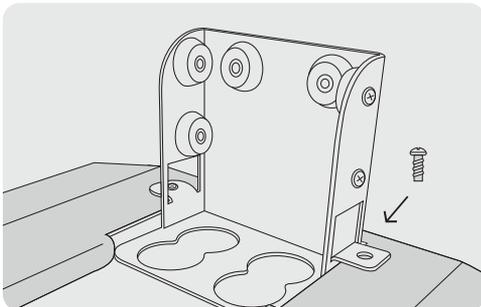


Step 4:

Stand assembly complete.

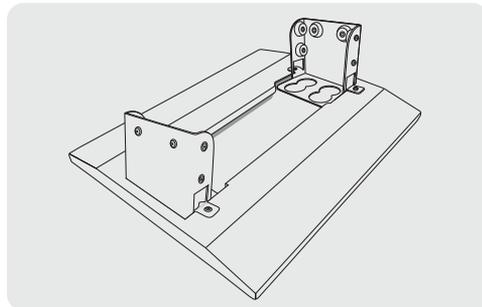

Step 5:
Completed look.

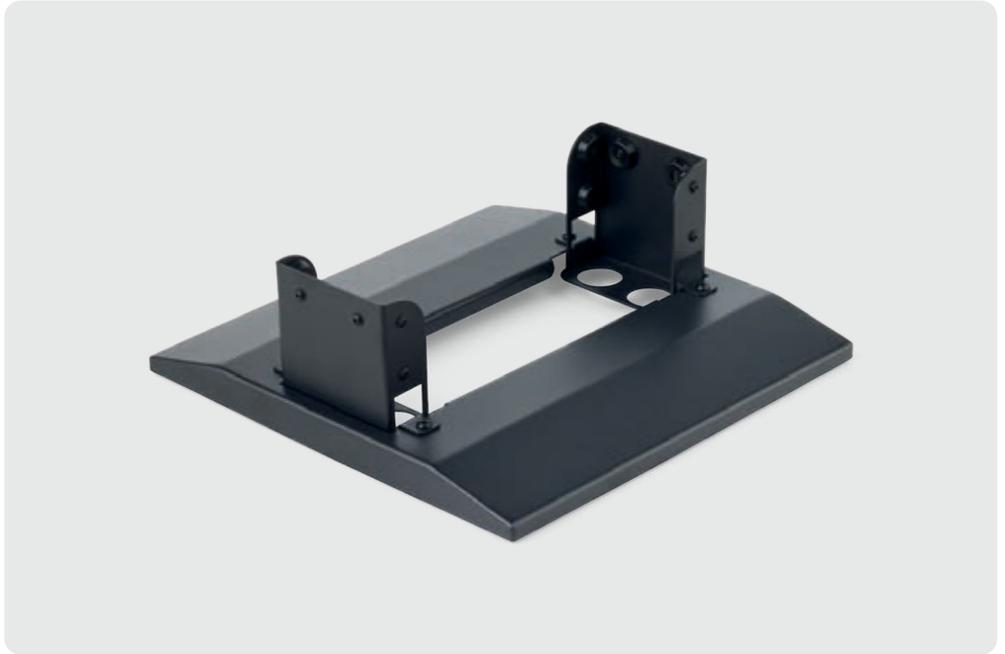

Step 6:
Place your Max or Super Max Red Light Therapy Device into the stand until it's firmly in place.





BON CHARGE

 @boncharge

 @boncharge_

 boncharge.com